



# TGA Balance Purge Switch

For Simple Gas Changeover for TGA Without  
Using Any Tools or Undoing Connections

[www.lab-ware.co.uk](http://www.lab-ware.co.uk)



Changing gases to any analyser usually involves taking out some tools and opening connections to change to a different gas. This can involve shutting down the gas system to avoid wasting gas, or may cause leaks when reconnected. The RSL gas switch removes the need to change any connections, and instead a simple changeover valve switches the gas to the new configuration. The existing second gas outlet is not affected by the switch, although different configurations can be provided as a customisation at addition cost.

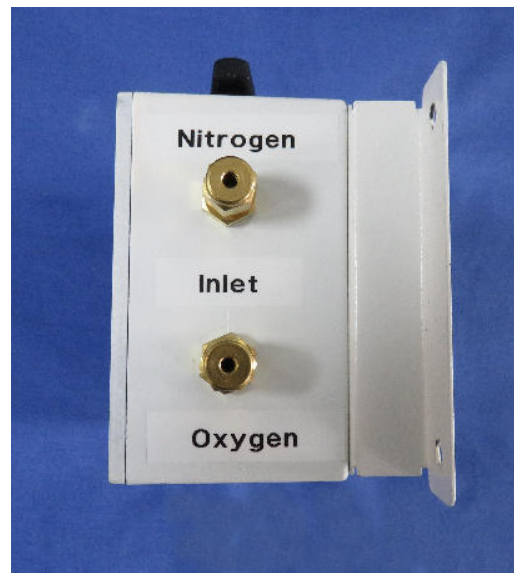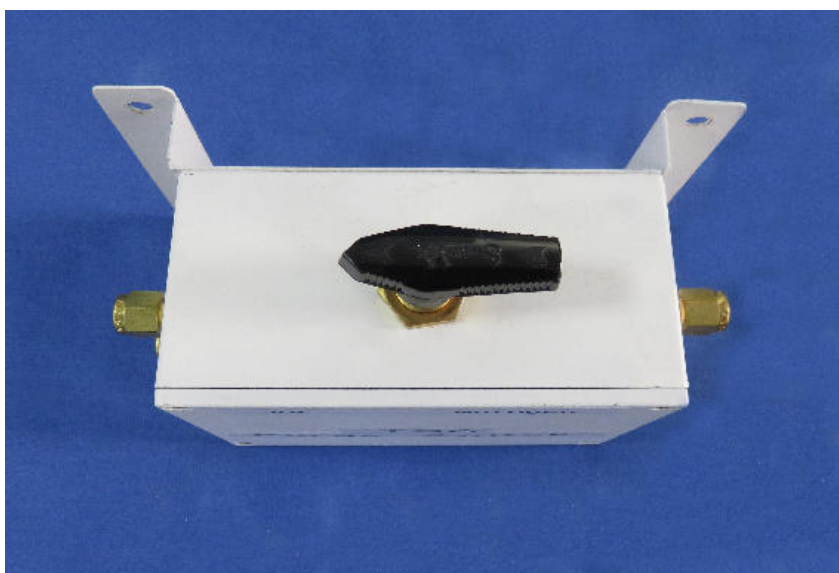

# TGA Balance Purge Switch

www.lab-ware.co.uk

## Specifications

|                       |  |
|-----------------------|--|
| Case                  | All-Metal case with four 5mm wall mounting holes at 76mm x 144mm centres.                    |
| Inlet Connections     | Brass Swagelok™ connection for 1/8 inch pipe connections                                     |
| Outlet Connections    | Brass Swagelok™ connection for 1/8 inch pipe connections                                     |
| Internal Pipework     | All in metal tubing to avoid gases diffusing from outside the pipes                          |
| Dimensions            | 120mm x 165mm x 78mm (W x H x D)   |
| Shipping Weight       | 0.70 kg  |
| Configuration Options | Different configurations of the inlet and outlet connections are available - please inquire. |

## Ordering Information

|                    |   |
|--------------------|---|
| Quotation Requests | Please email your requirements to <a href="mailto:sales@rustons.org.uk">sales@rustons.org.uk</a> , stating what you would like to order<br>Quotations are normally valid for 60 days.   |
| Ordering           | Please send your official order number to the address below or by email to <a href="mailto:sales@rustons.org.uk">sales@rustons.org.uk</a><br>If your organisation is VAT registered, please ensure that your VAT number is clearly visible on your order. |
| Delivery           | Delivery is normally 4 - 6 weeks from receipt of order. We may be able to provide faster delivery in some circumstances   |

